



# Rejeunesse™

Innovation of Dermal Rejuvenation



# Rejeunesse™

Innovation of Dermal Rejuvenation

Rejeunesse™ with lidocaine is a sterile, transparent, pyrogen-free, viscoelastic gel of 100% crosslinked hyaluronic acid of non-animal origin containing 0.3% lidocaine hydrochloride for treatment of wrinkles and folds. The lidocaine provides the patient's pain reducing effect during treatment.



Rejeunesse™ Sparkle

Rejeunesse™ Fine

Rejeunesse™ Deep

Rejeunesse™ Shape

Rejeunesse™ Contour




## 3 values of Rejeunesse™

Non-animal hyaluronic acid & no residual crosslinker



Long-lasting

More stable particle forms of HA, 100% crosslinked HA molecules

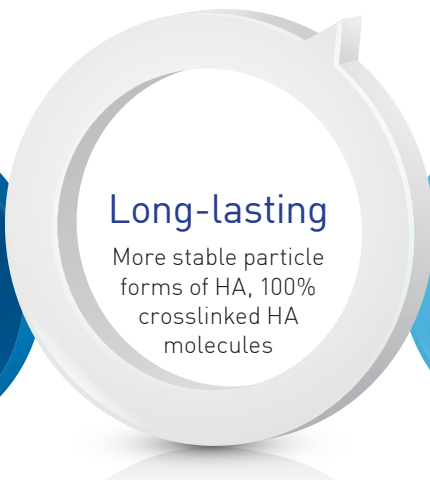

5 type molecular size & various applications



# Kind of Products



## Rejeunesse™ Sparkle

Rejeunesse™ Sparkle with lidocaine is crosslinked dermal filler for skin rejuvenation via superficial dermis injection.

## Rejeunesse™ Fine

Rejeunesse™ Fine with lidocaine is used in the correction of fine to medium sized depression of the face via mid to superficial dermis injection.



## Rejeunesse™ Deep

Rejeunesse™ Deep with lidocaine is used in the correction of medium to deep sized depression of the face via mid to deep dermis injection.

## Rejeunesse™ Shape

Rejeunesse™ Shape with lidocaine is used in the correction of deep sized depression of the face via deep dermis to subcutaneous injection.



## Rejeunesse™ Contour

Rejeunesse™ Contour with lidocaine is specially used in facial contouring, chin, nose end and shape augmentation of the face via deep dermis to subcutaneous injection.

# Rejeunesse™

Innovation of Dermal Rejuvenation




## Level of Injection

Rejeunesse™ Sparkle  
Rejeunesse™ Fine

Rejeunesse™ Deep

Rejeunesse™ Shape  
Rejeunesse™ Contour

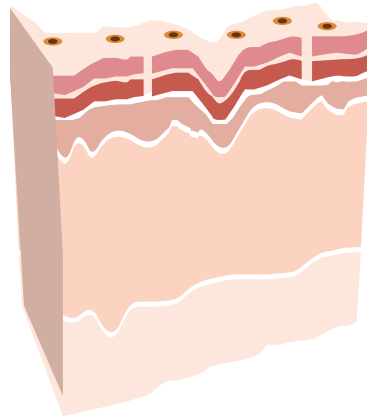

Superficial to Top Dermis  
Filling

Mid to Deep Dermis/Subcutaneous Layer  
Filling, Augmentation, and Contouring

Deep Subcutaneous Layer  
Filling, Augmentation, and Contouring



## Complete Product Line

	Rejeunesse™ Sparkle (with Lidocaine)	Rejeunesse™ Fine (with Lidocaine)	Rejeunesse™ Deep (with Lidocaine)	Rejeunesse™ Shape (with Lidocaine)	Rejeunesse™ Contour (with Lidocaine)
HA Concentration	20 mg/mL	24 mg/mL	24 mg/mL	24 mg/mL	24 mg/mL
Lidocaine	0.3%	0.3%	0.3%	0.3%	0.3%
Indications	Hydration and rejuvenation of skin	Superficial wrinkles, Crow's feet,	Deep skin folds, facial contouring	Facial contouring, chin and cheek augmentation	Facial contouring, chin, nose end & shape augmentation
Injection Site	Superficial dermis	Superficial to mid dermis	Mid to deep dermis	Deep to subcutaneous	Deep to subcutaneous
Packaging	3 x 2.5 mL	1 x 1.1 mL	1 x 1.1 mL	1 x 1.1 mL	1 x 1.1 mL
Needle Size	-	30G	27G	26G	23G
Duration	6-8 months	8-12 months	12-18 months	12-18 months	16-20 months
Storage	2-25°C	2-25°C	2-25°C	2-25°C	2-25°C
Shelf Life	24 months	24 months	24 months	24 months	24 months

# Rejeunesse™

Innovation of Dermal Rejuvenation

## Safety

Rejeunesse™ is non-animal based, crosslinked dermal filler using BDDE for linking agent. Rejeunesse™ product is obtained through our unique UPHEC technology, making it safe and non carcinogenic for humans. Rejeunesse™ contains less than 1.2EU/mL of endotoxins with no BDDE residue detected and it exceeds the highest level of HA purity. Rejeunesse™ is completely free of all animal products thus reducing the risk of an immunogenic reaction.

## Long-lasting

Once injected, HA is gradually broken and absorbed into the body, which is why how strongly HA is crosslinked is crucial to the longevity of the dermal fillers. With unique crosslinking technology, Rejeunesse™ has longevity over 12months.

## Efficacy

Rejeunesse™ offers visible results immediately and stable volume effect with minimized swelling after application.

## Treatment without pain

Rejeunesse™ has optimized the product to minimize the pain during and/or post injection for patients. Monophasic Rejeunesse™ enables smoothest and stable injection during procedure by maintaining consistent injection force. Especially, Rejeunesse™ Lidocaine series which combines 0.3% of Lidocaine (local anesthetic) minimizes injection pain during the procedure. Thus, it allows patients to experience the most comfortable, almost painless, procedure of dermal fillers.



Rejeunesse™



韩国海昇(株)  
KOREA HAESEUNG CO.,LTD.